

The Comfort hose manifold – an individual solution to meet all requirements



COMPRESSED AIR
SUCTION
IMAGING
DENTAL CARE
HYGIENE

More variation for more opportunities

The Comfort hose manifold by Dürr Dental adapts to every single working method

The equipment used with hose manifolds varies from surgery to surgery depending on the main treatments being carried out, the dentist's preferred suction technology, and handling methods. The Dürr Dental Comfort hose manifold has a modular design that allows it to be expanded to meet individual requirements. It is designed for up to five different components.

Whether you use a large or small suction hose, air/water spray, intraoral camera, speculum, or filter unit, you will find a place for it here. The vacuum-controlled selective manifold always provides the hose being used with the required suction power – with no disturbing noises. The Dürr Dental Comfort hose manifold can be mounted directly on the wall or – optionally – can be attached below the work surface with a pull-out pivot slide.

Camera bracket
Large suction hose/
spray mist suction system

Air/water spray

Small suction hose/
saliva ejector

Pivot slide



Really convenient filter replacement.

The filter element has a high collection volume, which means that the disposable filters have a long lifetime. When the filter does need to be changed, it's a convenient and hygienic process: Simply pull out the withdrawal rod along with the filter by the handle and throw the whole lot away.

This is what counts:

- Modular design, can be enhanced individually and easily
- Selection control for large and small suction hoses
- Pivoting hose connections on the filter housing to ensure a large range and low stress on the suction hoses
- Extremely ergonomic thanks to hose holder tilt adjustment function
- Optional shortening set for individual modification of hose length
- Can be mounted directly on the wall or, optionally, attached beneath the work surface with a pull-out pivot slide
- Long filter lifetime

Custom handpieces

The Dürr Dental handpiece range

Dürr Dental has the right hose adapter for practically every single manufacturer on the market. This enables the made-to-measure configuration of the respective suction handpieces. Users who want to work ergonomically and without getting tired should use the swivel joint. It ensures that the angle is always correct, reduces unnecessary leverage, and significantly improves handling. Each handpiece is available with and without a rotary slide. The slide can be used to temporarily interrupt the suction power. The rotating cover with Protect auxiliary air inlets further facilitates work.


This is what counts:

- Tailored configuration of suction handpieces
- Rotary slide for interrupting the suction power
- No disturbing flow noises when suction hoses are mounted, thanks to the use of special valves
- Swivel joint for high operating comfort (optional)
- Handpieces can be sterilized in an autoclave at up to 134 °C
- Swivel joint increases the lifetime of the suction hose

Example: Combination 1



Example: Combination 2

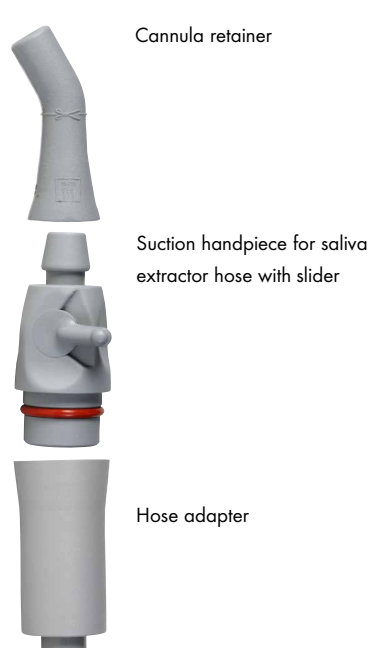

Example: Combination 3



Our range of suction cannulas



- Grey Pack of 20 0700-059-00
- Grey Pack of 5 0700-059-50
- Yellow Pack of 5 0700-059-51
- Blue Pack of 5 0700-059-54

Against reflux: Protect universal cannula

Secondary air inlets can prevent the reflux of contaminated secretion from suction system when suction is applied to the cannula. The Protect universal cannula has no corners or edges and leaves no sore spots.



- Grey Pack of 4 0700-058-50

Prophylaxis cannula: with up to 2x more aerosol uptake

Only at Dürr Dental: a specially developed rotatable and lockable protective shield enables increased aerosol uptake and particularly effective, ergonomic suction as part of prophylaxis treatments with powder jet cleaning devices.



- Grey Pack of 20 0700-054-00
- Grey Pack of 5 0700-055-50
- Yellow Pack of 5 0700-055-51
- Pink Pack of 5 0700-055-53
- Blue Pack of 5 0700-055-54
- Turquoise Pack of 5 0700-055-55



For adults: universal cannula

Dürr Dental universal cannulas are particularly quiet. In a study by the Fraunhofer Institute, the measured volume was significantly below the level of comparable competitor products. They therefore ensure certified, low-noise suction of spray mist. And all of this without leaving behind any pressure sores.




- Grey* Pack of 5 A700056050
- Yellow* Pack of 5 A700056051
- Pink Pack of 5 A700056053
- Blue* Pack of 5 A700056054
- Turquoise Pack of 5 A700056055
- ∅ 16 mm

For children: Petito universal cannula

Children require a highly functional cannula that permits ergonomic treatment. The Petito universal cannula is smaller than the universal cannula; its Protect secondary air inlets provide optimal protection against reflux from the suction system. It has no corners or edges and leaves no sore spots.

*Also available in ∅ 11 mm without Protect secondary air inlets

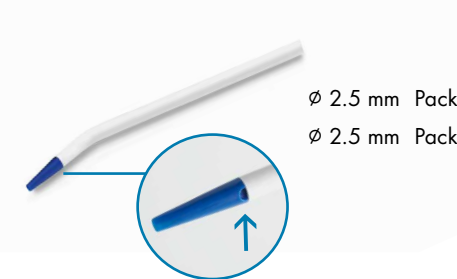

- ∅ 2.5 mm Pack of 20 0700-007-50
- ∅ 2.5 mm Pack of 100 0700-007-51

Surgical suction cannula for single use

Surgical cannula for single use with a 2.5 mm diameter and a cannula length of approx. 17 cm. The cannula has a tip with integrated Protect secondary air inlets. This offers the best possible protection from reflux from the suction system. The cannulas are individually packaged and sterile.

DÜRR DENTAL SE
 Höpfigheimer Str. 17
 74321 Bietigheim-Bissingen
 Germany
www.duerrdental.com
info@duerrdental.com

**DÜRR
 DENTAL**
 THE BEST, BY DESIGN